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Saljol Adapter Set for Seat Height Increase and Seat Angle Adjustment for Shower Stool & Shower Chair

Instruction manual

Manufacturer: Saljol GmbH
Product: Saljol Adapter Set for Seat Height Increase and Seat Angle Adjustment
Type: Accessories
Article number: 4058789017050

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Scope of delivery

4 adapter tubes, height 4 cm	1 Instruction manual
4 adapter tubes, height 8 cm	
4 M6 socket head cap screws	Required tools:
4 M6 nuts	5 mm Hex key
4 curved washers	10 mm open-ended spanner

Restrictions on use

The maximum permitted user weight is 150 kg.


Accident prevention

- Be aware of potential pinch points during assembly.
- The adapter tubes may only be fitted to the designated mounting points on the shower stool or shower chair.
- The adapter set may only be used with Saljol shower stools and shower chairs.
- Only use the adapter tubes, screws and washers supplied.
- Tighten all screws evenly and securely.
- Before each use, check that the adapter tubes are correctly positioned and that all screws are securely tightened.

Assembly

1. Remove the seat

- Slide the two seat sections one after the other evenly inwards in the direction of the arrow until they disengage from the mounts.
- Remove the seat sections completely from the frame.


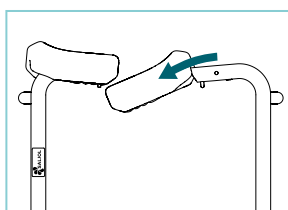
2. Fit the adapter tubes

First select the required seat height or seat angle adjustment and seat width:

- Increase seat height by 4 cm: Fit four adapter tubes, each 4 cm long.
- Increase seat height by 8 cm: Fit four adapter tubes, each 8 cm long.
- Seat angle of 13°: Fit either two 4 cm adapter tubes to the rear mounting points or two 4 cm adapter tubes at the front and two 8 cm adapter tubes at the rear (the side with the open frame is the front; the side with the closed frame is the rear).
- Seat width 36 cm: Fit four adapter tubes to the four designated mounting points **1**.
- Seat width 39 cm: Fit four adapter tubes to the four designated mounting points **2**.
- Seat width 42 cm: Fit four adapter tubes to the four designated mounting points **3**.


- Position the selected adapter tubes at the designated mounting points 1/2/3 on the frame.
- Align the holes in the adapter tubes with the holes in the frame.
- Insert the supplied M6 socket head cap screws together with the curved washers through the holes.
- Fit the nuts onto the screws.
- Tighten all fitted screw connections evenly using suitable tools.
- Then check that the adapter tubes are securely fitted.

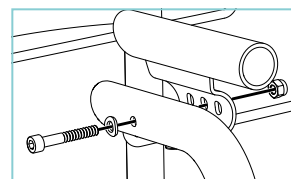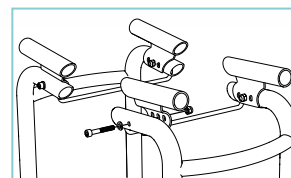3. Fit the seat

- Position the seat sections at an angle above the mounts on the frame.
- Slide the seat sections outwards onto the frame in the direction of the arrow.
- Press the seat sections fully into the mounts.
- Then check that the seat is securely and firmly fitted.
- Before use, ensure that all screw connections are securely tightened and that the seat of the shower stool or shower chair is secure and stable.

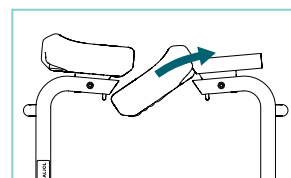Note

Tighten the screws only hand-tight to avoid damaging the material.

Technical Data

Seat height adjustment	54 cm/58 cm
Seat angle adjustment	13°
Hygiene opening, setting 1/2/3	6/9/12 cm
Seat width 1/2/3	36/39/42 cm

Cleaning and Disinfection

Important

- Do not use aggressive cleaning agents such as solvents or abrasive cleaners, or hard brushes.
- Do not clean the shower stool/shower chair using a high-pressure cleaner or steam cleaner.
- When disinfecting, follow the instructions for use and application provided by the manufacturer of the cleaning and disinfecting agents.
- After cleaning, rinse the shower stool/shower chair with clean, warm water and dry it thoroughly.
- After use, thoroughly remove residues of personal care products such as liquid soap or shower gel with clean, warm water

Formalities of the Instruction Manual

Version/revision: 2026.01

Creation date: 23/09/2026

Information was current at the time of going to press. We reserve the right to make technical changes. Technical drawings and photos may differ from reality.

©2026 Saljol GmbH: all rights reserved. Reprinting, including extracts, is only permitted with the authorisation of Saljol GmbH.