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Saljol Backrest for Shower Stool

Instruction manual

Manufacturer: Saljol GmbH
Product: Backrest
Type: Accessories
Article number: 4058789017043

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Figure:
Shower stool
with backrest

Scope of delivery

- | | |
|-------------------------|----------------------|
| 1 backrest | 4 spacer sleeves |
| 2 EPP backrest elements | 1 instruction manual |
| 4 M6 pan head screws | 1 5 mm hex key |
| 4 curved washers | |

Restrictions on use

The maximum permitted user weight is 150 kg.


Accident prevention

- Be aware of potential pinch points during assembly.
- An incorrectly fitted or loose backrest may cause loss of balance and falls. Before each use, check that the backrest is securely fitted and in proper condition. If necessary, retighten the screws and correctly adjust the EPP backrest elements.
- The backrest may only be fitted to the designated fixing points on the shower stool.
- Use only the supplied screws, washers and spacer sleeves.
- Tighten all screws evenly and securely.
- The backrest is intended solely to provide support while seated. It must not be used as an aid for standing up or as a support aid.
- Do not lean against the backrest suddenly or with excessive force.
- Exposure to external sources of heat or cold, such as direct sunlight or low ambient temperatures, may alter the surface temperature of the backrest. Check the surface temperature before use.

Assembly

- Position the backrest at the designated fixing points on the underside of the shower stool frame.
- Align the holes in the backrest with the holes in the frame.
- Insert the supplied screws together with the curved washers through the holes in the backrest, then place the spacer sleeves onto the screws.

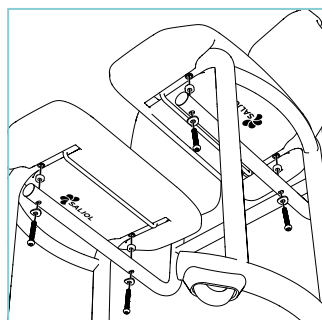
- Using the supplied hex key, tighten the screws evenly until hand-tight.
- Then check that the backrest is securely and firmly fitted.

Note

Tighten the screws by hand only to prevent damage to the material.

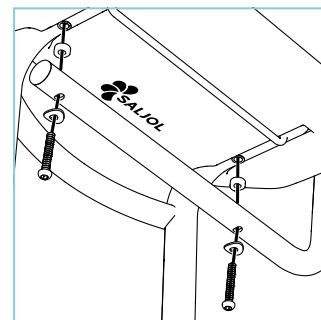Technical data

Recommended user height	150–190 cm
Maximum user weight	150 kg
Total weight	1,5 kg
Backrest length	37 cm
Backrest width	33/36/39 cm

Cleaning and disinfection

Important

- Do not use aggressive cleaning agents, such as solvents or abrasive cleaners, and do not use hard brushes.
- Do not clean the backrest or shower stool using a high-pressure cleaner or steam cleaner.
- When disinfecting, observe the instructions for use and processing provided by the manufacturer of the cleaning and disinfecting agents.
- After cleaning, rinse the shower stool with clean, warm water and dry it thoroughly to extend its service life.
- After use, thoroughly remove residues of personal care products, such as liquid soap, shampoo or shower gel, using clean, warm water.

Formalities of the Instruction Manual

Version/revision: 2026.01

Creation date: 28/08/2026

Information was current at the time of going to press. We reserve the right to make technical changes. Technical drawings and photos may differ from reality.

©2026 Saljol GmbH: all rights reserved. Reprinting, including extracts, is only permitted with the authorisation of Saljol GmbH.

EN