# Saljol Schwenktisch

Bedienungsanleitung | Instruction manual | Manuel d'utilisation  
Istruzioni per l'uso | Gebruiksaanwijzing | Manual de instrucciones

DE

EN

FR

IT

NL

ES

Vorwort

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Saljol  
Swivel Table

Foreword

Thank you for choosing the Saljol swivel table! With your new swivel table, you will enjoy greater comfort, flexibility and independence in everyday life. The stable and smooth swivel mechanism has been developed with a particular focus on safety and user-friendliness. Please read this instruction manual carefully before assembling and using the swivel table for the first time. It contains important safety instructions and practical tips for use and care.

EN

Contents

|    |                                                    |    |
|----|----------------------------------------------------|----|
| 1. | Identification .....                               | 7  |
| 2. | Scope of delivery .....                            | 8  |
| 3. | Notes.....                                         | 8  |
|    | 3.1 Safety instructions.....                       | 8  |
|    | 3.2 Instructions for safe use of the product ..... | 8  |
| 4. | Assembly instructions .....                        | 8  |
|    | 4.1 Preparation .....                              | 8  |
|    | 4.2 Mounting the fixing plate .....                | 9  |
|    | 4.3 Attaching the tabletop .....                   | 9  |
|    | 4.4 Height adjustment.....                         | 10 |
| 5. | Cleaning and disinfection.....                     | 10 |
| 6. | Storage and transport.....                         | 10 |
| 7. | Technical data .....                               | 10 |
| 8. | Warranty .....                                     | 11 |
| 9. | Disposal .....                                     | 11 |

1. Identification

Manufacturer: Saljol GmbH  
Product: Swivel table  
Type: Side table  
Version: 2025  
Article number: 4058789012864

Manufacturer:  
Saljol GmbH  
Frauenstraße 32, 82216 Maisach, Germany  
Phone: +49 8141 317740, Fax: +49 8141 3177411  
Email: [hallo@saljol.de](mailto:hallo@saljol.de), Web: [www.saljol.de](http://www.saljol.de)

Swiss authorised representative and importer:

|    |     |                  |
|----|-----|------------------|
| CH | REP | CON GmbH         |
|    |     | Talgut-Zentrum 7 |
|    |     | CH-3063 Ittigen  |

Information on this instruction manual  
Version/Revision: 2026.01  
Date of issue: 04/02/2026

The information corresponds to the status at the time of printing. Technical changes reserved. Technical drawings and photographs may differ from the actual product.  
©2025 Saljol GmbH. All rights reserved. Reproduction, even in part, is permitted only with the consent of Saljol GmbH.

## 2. Scope of delivery

- 1 tabletop
- 1 felt mat
- 1 support tube with mounted fixing plate
- 3 cross-head thumb screws M5×11
- 1 hex nut M8
- 2 countersunk Allen screws M10×30
- 2 spacers for height adjustment for M10×30 screws
- 1 instruction manual

## 3. Notes

### 3.1 Safety instructions

This instruction manual uses safety instructions that are identified in two ways: by a symbol and by a signal word. The symbol indicates the type of instruction. The signal word provides further explanation, for example the severity of a potential risk.

#### Caution / Important

Potentially hazardous situation (minor injuries or material damage), obligation to observe specific behaviour or actions for safe handling of the swivel table.

#### Note

Application tips and particularly useful information.

### 3.2 Instructions for safe use of the product

For safe use of the product, please observe the following:

#### Restrictions on use

- The support tube can bear a max. load of 100 kg.
- The tabletop may be loaded with a max. of 20 kg.
- Never use in the presence of flammable or burning substances.
- Do not lean on, sit on or climb onto the tabletop.
- Allow children near the swivel table only under supervision.
- The swivel table is intended for indoor use only.

#### Caution

- Keep the swivel table free from obstacles to avoid trapping or impact.
- Do not use the swivel table in damp or wet environments.
- Do not place hot, liquid or unstable objects on the tabletop.
- The swivel table may only be attached to suitable riser recliner chairs and only to the lower frame.


## 4. Assembly instructions

### 4.1 Preparation

- You will need a cross-head screwdriver and a 6 mm Allen key.
- First disconnect the chair from the power supply.
- It is recommended to assemble the table on the chair with two people.
- Ensure that no cables, objects or animals are near the chair.
- The side table can be mounted on the left or right side of the chair. Decide which side it should be mounted on.
- Lay the chair on the side on which the table will not be mounted. To do so, hold the chair by the headrest and side parts and carefully turn it onto its side.


In this description, the table is mounted on the left when viewed from the front. Therefore, the chair is laid on its right side.

- Now unscrew the marked foot from the frame.

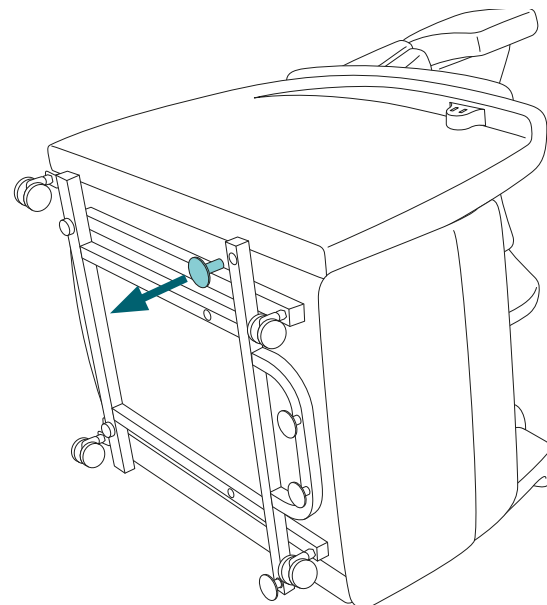
### 4.2 Mounting the fixing plate

- Place the fixing plate against the chair frame. The support tube may rest on the chair.
- Prepare the two countersunk Allen screws M10×30 and the two spacers for height adjustment.
- Lift the support tube with one hand so that the fixing plate can move freely.
- Place the spacers between the fixing plate and the chair frame at two diagonally opposite holes (see illustration).
- Tighten the two screws alternately using the Allen key. This ensures even fastening of the fixing plate and prevents jamming.

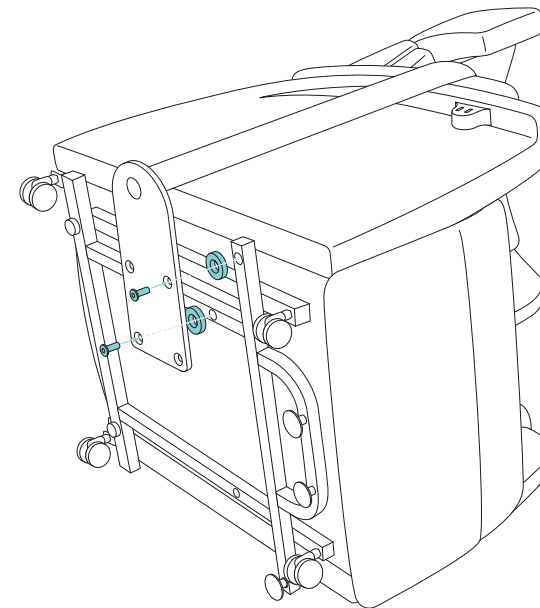
- Then reattach the previously removed foot using the M8 hex nut.

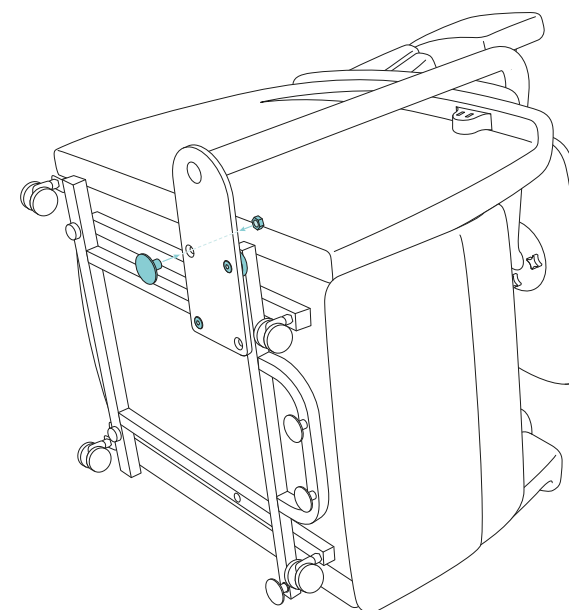- All parts marked in the illustration should now be firmly tightened.

### 4.3 Attaching the tabletop

- Take the tabletop and prepare the three cross-head thumb screws M5×11.
- Place the tabletop onto the support tube – the drill holes in the tabletop should be visible through the long slots of the mounting plate.
- Initially screw the three thumb screws loosely into the visible holes.
- Align the tabletop to the desired position.
- Then tighten all screws firmly.
- Now return the chair to the upright position.

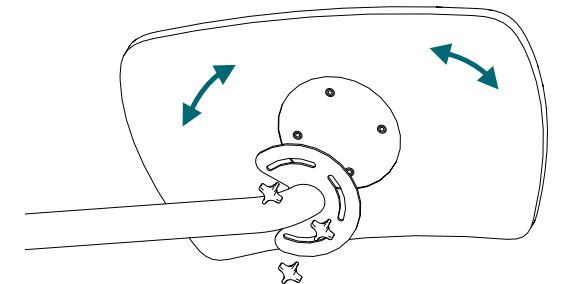
#### Note

The maximum load at the outer edges of the tabletop is 20 kg. The maximum load on the support tube (e.g. when leaning on it) is 100 kg.

#### Important

Always swivel the table outwards before activating the stand-up function. There is a risk of crushing if the chair is raised without the swivel table being swung to the side. The support tube of the table can also be used as a standing aid. Hold onto the area marked in petrol.

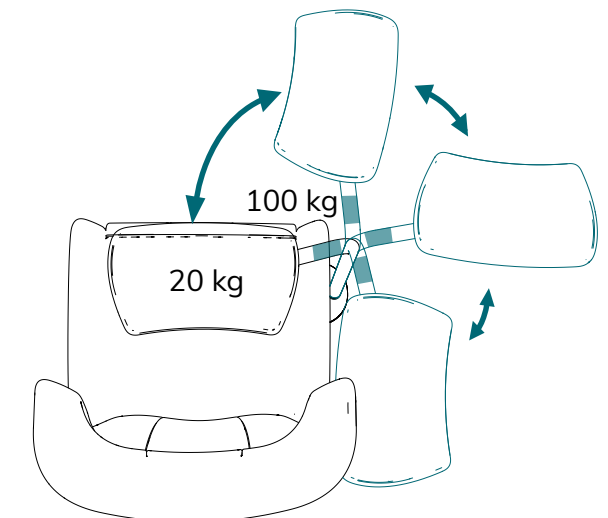
 **Important**

It is prohibited to lean on the tabletop, as this poses a risk of falling.

**4.4 Height adjustment**

- The swivel table can be adjusted to three heights: 70.5 cm, 73.5 cm and 76.5 cm. The table is supplied with a height of 76.5 cm.
- Ideally, work with two people when adjusting the height. Choose the height that offers the greatest comfort to the user.
  - Lift the support tube to gain access to the extension parts mounted on the base.
  - Loosen the two Torx countersunk screws of the upper extension part.
  - Remove the loosened extension part to reduce the table height by 3 cm.
  - If further reduction is required, repeat the procedure and remove the second extension part as well.
  - After removing the parts, the support tube automatically lowers towards the fixing plate.
  - Keep the removed extension parts and screws in case they need to be reinstalled.

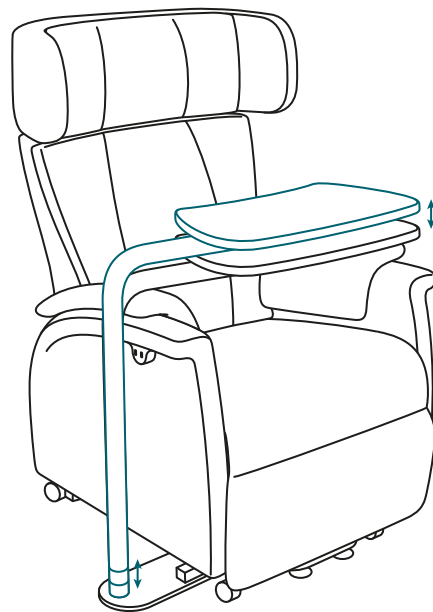**5. Cleaning and disinfection**

Clean the swivel table according to the degree of soiling, but at least once a month. Use a soft cloth or sponge and warm water. For stubborn dirt, a commercially available cleaning agent may be added. The table can be disinfected using household disinfectants. Always follow the application instructions of the cleaning and disinfecting agents used. To eliminate occasional squeaking or creaking of the swivel mechanism, simple products such as PTFE (Teflon) spray may be used.

 **Important**

Never use aggressive cleaning agents such as solvents, abrasive cleaners or hard brushes.

**6. Storage and transport**

Temperature: 0 °C to +35 °C  
Humidity: 45 % to 60 %

When not in use, store the table in a dry place inaccessible to children or unauthorised persons. If the table is to be stored or transported, it is advisable to keep the original packaging so that it is available if required.

**7. Technical data**

|                             |                                     |
|-----------------------------|-------------------------------------|
| Tabletop dimensions         | 52 x 30 x 1.5 cm                    |
| Height above floor          | 70.5 to 76.5 cm                     |
| Maximum load (table edge)   | 20 kg                               |
| Maximum load (support tube) | 100 kg<br>(e.g. when leaning on it) |
| Tabletop                    | Birch wood                          |
| Tube                        | Steel, KTL-coated                   |
| Felt mat (not washable)     | 90 % wool + 10 % polyester          |

**8. Warranty**

If the purchased item is defective, the statutory warranty provisions apply. If you purchased the product from a specialist retailer, please contact your seller immediately. If purchased directly from Saljol, please contact Saljol directly at +49 8141 317740. Obvious defects must be reported in writing within ten days of delivery, hidden defects within ten days of discovery. These periods are exclusion periods. Saljol is entitled to decide whether a defect will be remedied by repair or replacement. In all other respects, the General Terms and Conditions of Saljol apply.

**9. Disposal**

Dispose of the materials separately by type in the designated collection systems. Dispose of hazardous waste via the appropriate collection systems.

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