



# WR

Indoor Rollator



The sturdiest Indoor Rollator  
in the world – safely practical

## The Saljol Indoor Rollator

### Premium quality with many extras

This premium product has been specially developed for your safety at home. You can rest on the seat board or move through the rooms. The directionally lockable rear wheels make it easier to drive straight ahead. You can use the brake bar, which matches the colour of the rollator frame, as an alternative to the brake handles. You can customise your rollator even further with additional accessories.

### Stable like a mobile handrail – optimal support with every movement

The indoor rollator provides extra support within your home, especially in cases of unsteady walking, dizziness, or balance issues. Thanks to its sturdy, ergonomically designed frame, you can safely lean on any side of the rollator at all times without the risk of it toppling over.








### Internal brake cables – innovative solution for your safety

All the brake cables of the rollator run through the sturdy frame, eliminating the risk of snagging. This ensures you can move safely through rooms without getting caught on furniture and place everyday items on the tray without getting entangled in the cables.



### Double the comfort when braking

Thanks to the brake bar, you can slow down the rollator evenly with one hand or activate the parking brake. You don't need to worry about where you grip the bar, as it can be operated from any position.

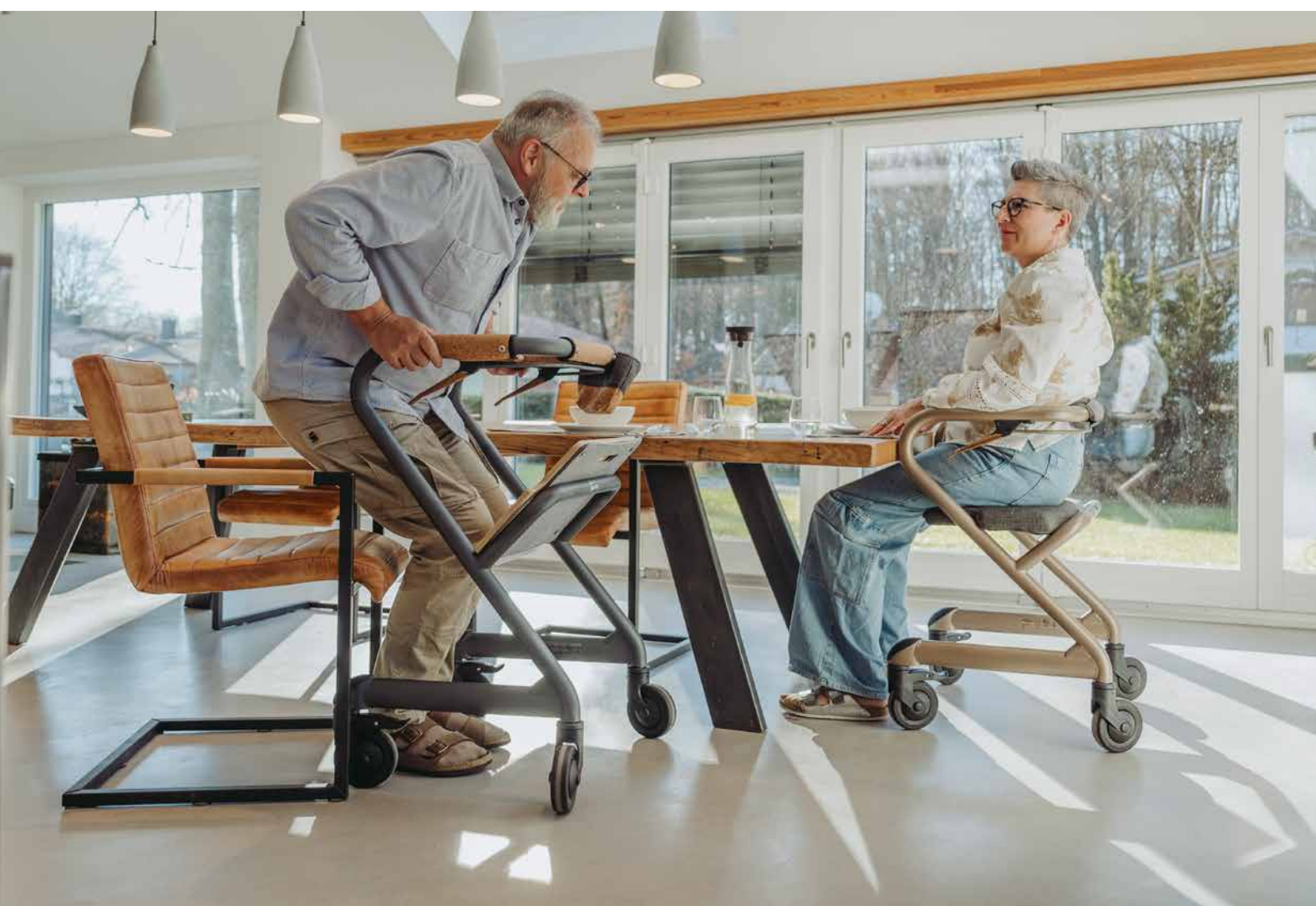

### Rolling seat – enables easy movement while sitting by tripping

In cases of rapid fatigue, fall risk, weak leg muscles, or joint instability, the rollator offers a large seat with a warm felt cushion. Tripping – moving forward while seated using your feet – is explicitly permitted with this rollator. Safely navigate through different rooms or comfortably prepare your meals in the kitchen while sitting.

### Lockable rear wheels – more control and stability when moving straight

The rear wheels can be locked in the longitudinal direction to provide more control and stability while moving straight. This feature is especially helpful when you need to lean more heavily on the rollator.








### Flexibility even in tight spaces

The indoor rollator ensures safe navigation through narrow doorways and cramped or cluttered hallways and rooms. With its 360° maneuverability, it can be easily moved sideways or lengthways through doors or tight spaces between furniture. Stability is maintained at all times, allowing you to lean on it securely from any direction.





### A highly practical bathroom aid: walk, stand, and turn on the spot

The indoor rollator offers reliable support for daily hygiene and bathroom visits. Thanks to its 360° maneuverability and thoughtfully designed frame, you can comfortably lean on it from the side or back and position it directly in front of the sink. It can easily roll over a standard-height toilet, allowing you to hold on securely while sitting down or standing up.

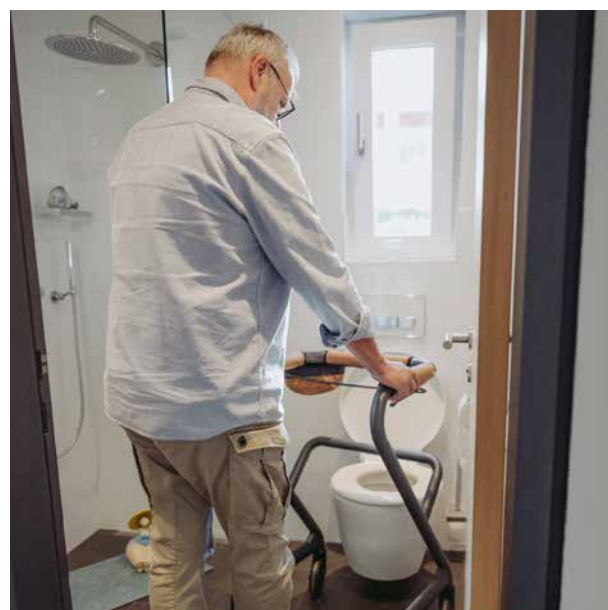

### Versatile household helper for every situation

The indoor rollator provides reliable support for transporting everyday items. Its sturdy tray with a non-slip felt pad offers a secure way to serve food and drinks, keeping cups and plates firmly in place. For larger items like laundry baskets or flower pots, the durable seat board serves as a practical surface. With your hands securely on the rollator, you can effortlessly and safely move everything from one place to another.

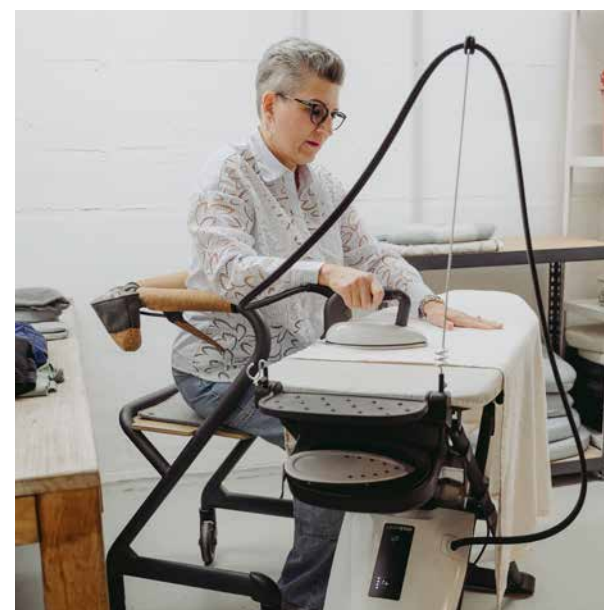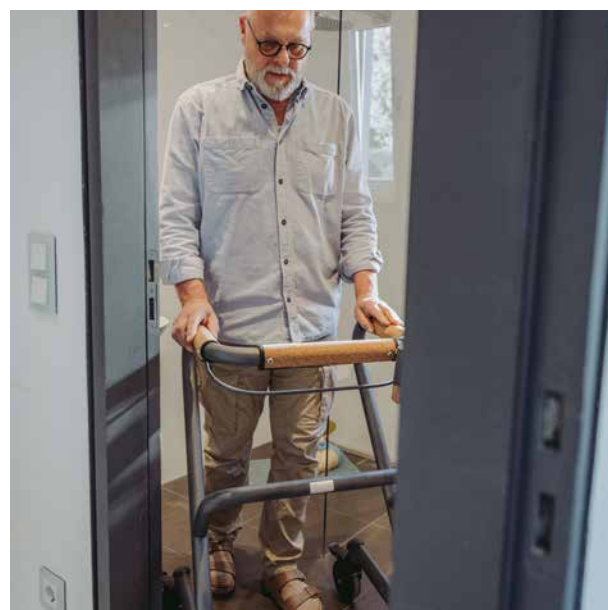





### **Practical bag – essentials always within reach**

The removable, colour-coordinated bag ensures that personal items like keys, a phone, or tissues are always close at hand. Whether you're walking without the tray – such as during a trip to the bathroom or while moving through your home – or using the tray for serving, the bag keeps your essentials easily accessible.




### **The right product for you – available for every need**

Saljol Indoor Rollators are available in four sizes, elegant colours, and additionally with a wider seat and higher weight capacity. This gives you the option to choose the solution that best suits your needs and personal taste. Whether for more comfort, functionality, or style – you will find the perfect combination to meet your requirements.



### **Skin-friendly handrail pads for maximum comfort**

The removable, skin-friendly handrail pads made from cork fabric and neoprene feel exceptionally comfortable. The shock-absorbing and heat-retaining materials support the mobility of your hands and fingers, especially if you have sensitive skin. Their non-slip, quick-drying properties ensure a secure grip, even in wet conditions.

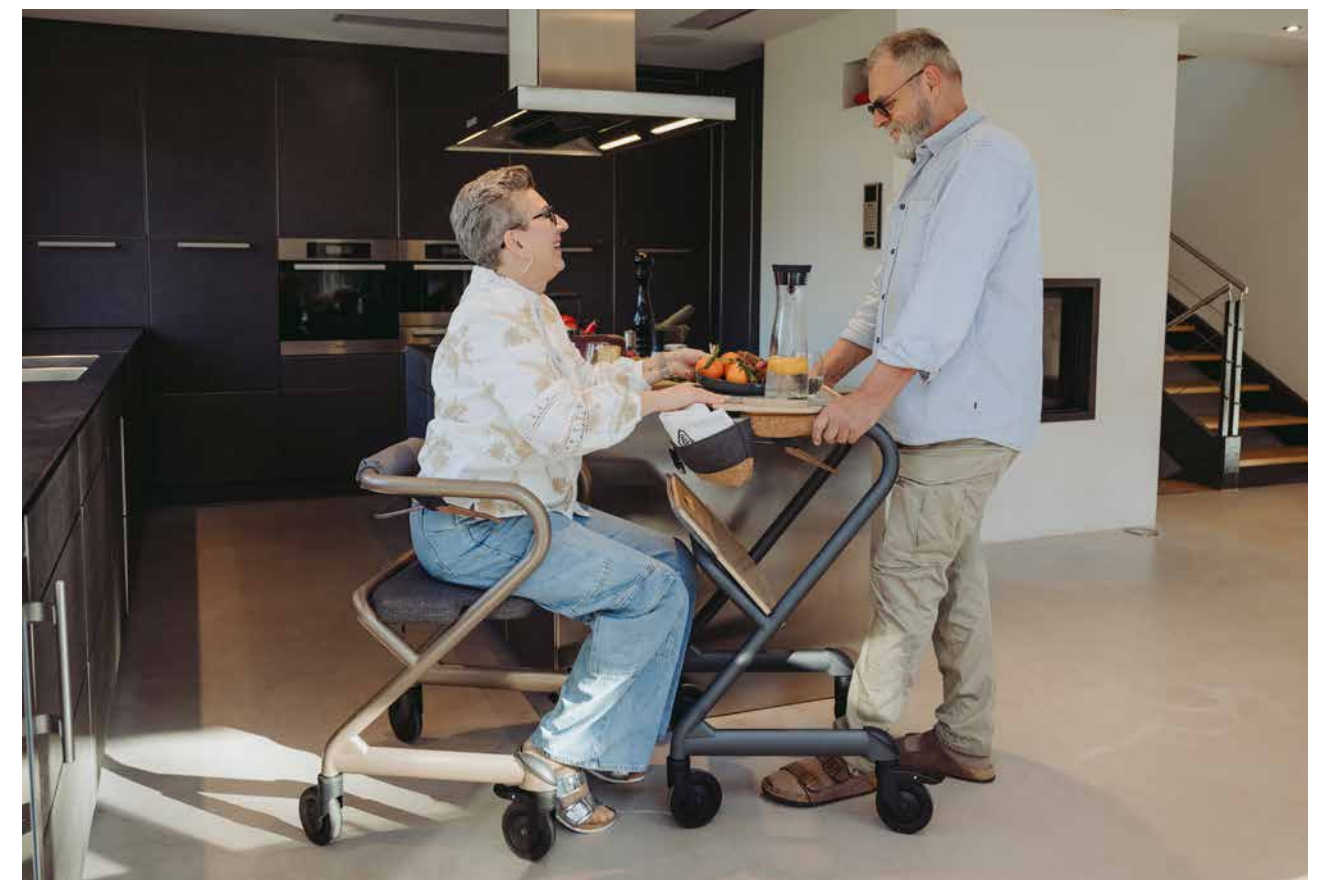





German Engineering from Bavaria

German Design



### Great luxury that remains affordable

Saljol Rollators offer the perfect mix of functionality, quality, and a touch of luxury. German engineering combines high-quality features with comfort to make your everyday life easier – all at a fair overall price.



### Comprehensive support – for your mobility and safety

Our physiotherapist, Katharina, supports you in all areas where expert physiotherapeutic knowledge is required. Our helpful customer service is always at your side – whether in the medical supply store or directly with us. If you have any questions or uncertainties regarding your Saljol Indoor Rollator, we provide the support you need to get the most out of your rollator.



### Training with your rollator

To avoid medical side effects when using our rollators, we invented Rollastic. It's a tailored exercise programme for rollator users that, when performed regularly, can help maintain mobility and strength. Special videos are provided with the rollator. For your safety in road traffic, we developed the safety training „Rollator Mobility“. Helpful tips and tricks can be found at [www.rollatormobility.com](http://www.rollatormobility.com).



### Claudius: “Being able to sit and move on the rollator makes my daily life so much easier!”

„Whether I’m cooking, unloading the dishwasher, or moving around my flat, the indoor rollator is a great support. Not as a typical rollator, but as a readily available seat. When I become exhausted and can’t stand anymore – which happens very quickly due to my MS – I can simply lean on it and sit down. That’s the biggest benefit for me. It’s like a mobile seat that allows me to move around my home with ease.”



Indoor Rollator	WR43	WR50	WR53	WR56	WRXL*
Recommended body height	135–155	145–165	150–180	170–200	170–210
Maximum user weight	120 kg	150 kg	150 kg	150 kg	180 kg
Handle height	78 cm	80 cm	86 cm	92 cm	95 cm
Seat height	43 cm	50 cm	53 cm	56 cm	58 cm
Width bottom (wheels)	56 cm	56 cm	56 cm	56 cm	68 cm
Width top	51 cm	51 cm	51 cm	51 cm	65 cm
Distance between handles	43 cm	43 cm	43 cm	43 cm	60 cm
Weight without accessories	9.4 kg	9.6 kg	9.8 kg	10.0 kg	13.4 kg
Compatible with accessories	✓	✓	✓	✓	✗

### Indoor Rollator accessories



Small bag



Tray



Back roll



Soft seat



JOST bag  
beige



JOST bag  
transparent

### Frame colour



Velvet Grey



Champagne Gold

\*WRXL is supplied with seat board and brake bar in frame colour and is only available in Velvet Grey.